






* Measure Length of Middle Finger

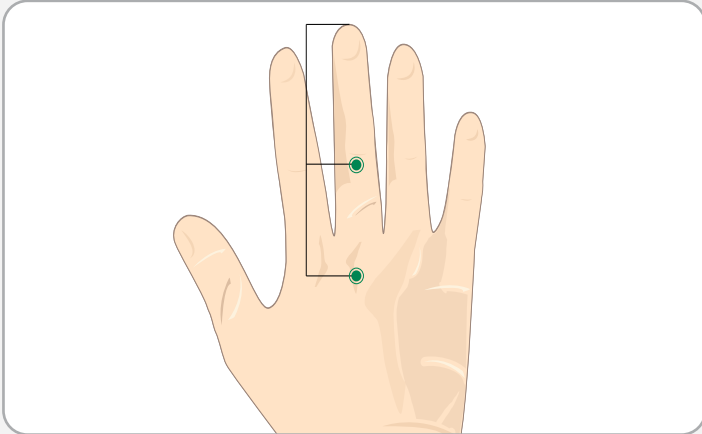


1 Locate Flexion Point & Calculate Distance

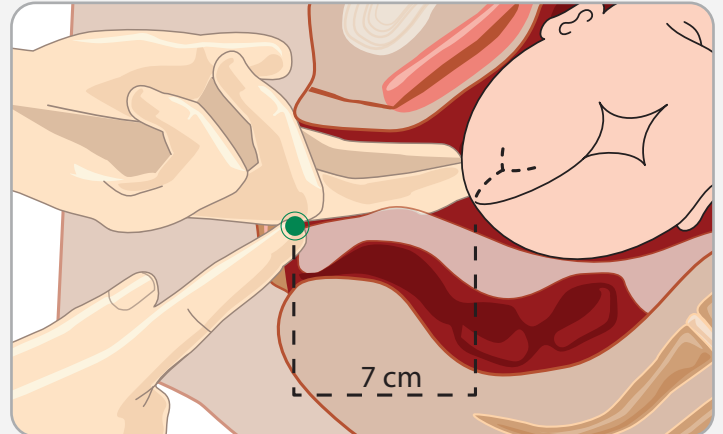


2 Hold & Insert Cup

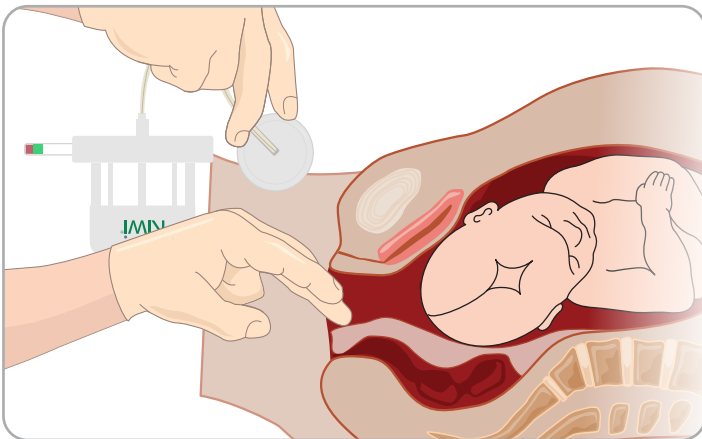


3 Move Cup Over Flexion Point

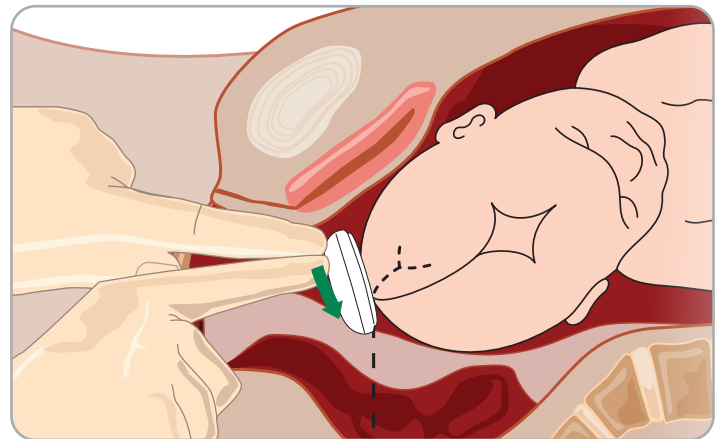

4 Create Vacuum & Exclude Maternal Tissue

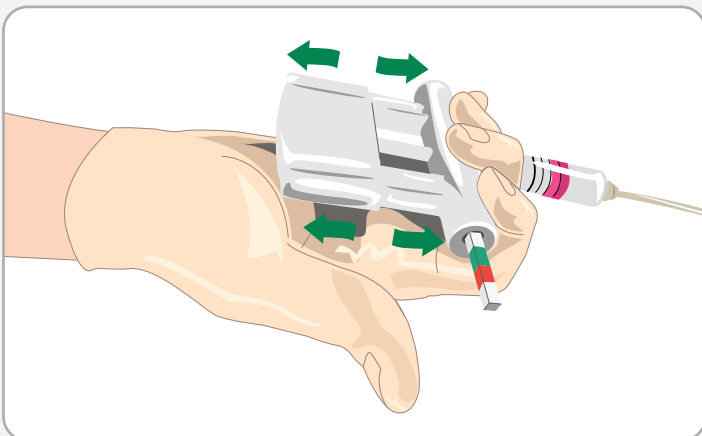

5 Initiate Traction

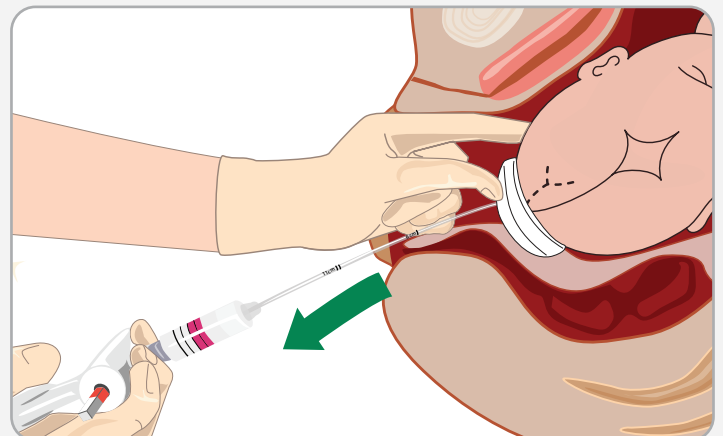


5 Step Key Takeaways

1 Locate Flexion Point & Calculate Distance

Locate the flexion point. With opposite hand, note distance where finger meets introitus.

This is the distance you will insert the cup in step 3.

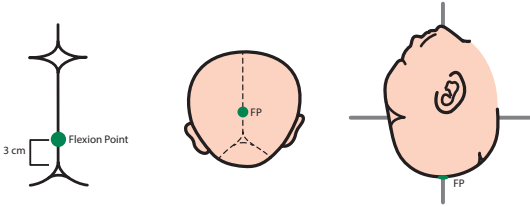


NOTE: The **flexion point** is located on the sagittal suture 3 cm forward of the posterior fontanelle.

2 Hold & Insert Cup

Hold cup with thumb on tube in groove, fingers on foam and insert into the vagina. After inserting, adjust cup so it lies flat against fetal head along midline.

NOTE: Place 'groove' and stem of cup at **12 o'clock**. This allows one to visualize rotation as descent occurs.

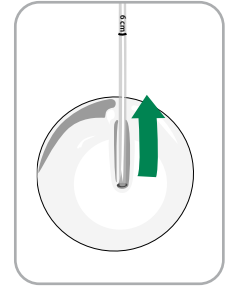


3 Move Cup Over Flexion Point

Note 6 & 11cm markings on stem tubing. Reference these markings to know how far to insert cup. Push cup **posteriorly along maternal midline** over flexion point. Insertion distance is the distance measured in Step 1.

NOTE: The flexion point is along the midline.


Do not push cup laterally.

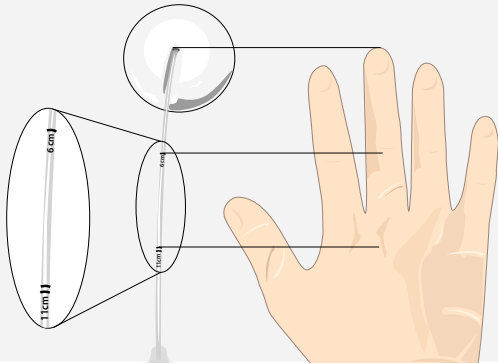

4 Create Vacuum & Exclude Maternal Tissue

Use palm pump to create vacuum. Once initial vacuum is created, feel around cup to exclude any maternal tissue. Then continue pumping vacuum to 600mmHg.

NOTE: Create vacuum to green section on gauge. If possible, 'bury' the green so it is no longer visible.

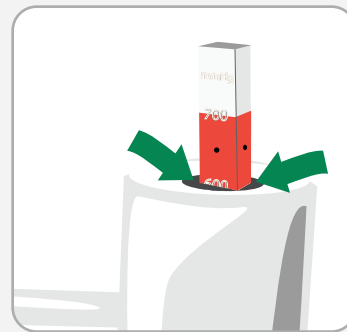


5 Initiate Traction

Initiate traction along axis of pelvis. Pull during contractions. **Do not move handle up and down or side to side while pulling.**

FINGER / THUMB TECHNIQUE


Use thumb and index finger of non-pulling hand to support the cup and the fetal head while pulling. This will help identify loss of vacuum before a pop-off occurs, while also providing feedback on descent.

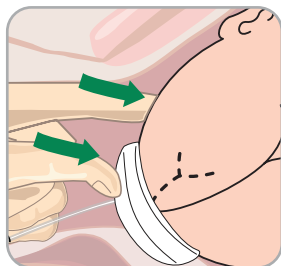


2 x 2 FINGERTIP GRIP

Ensure you have two fingers on the handle to each side of the stem. This helps with equal pulling, while reducing likelihood of accidentally pressing the vacuum release button.

Correct

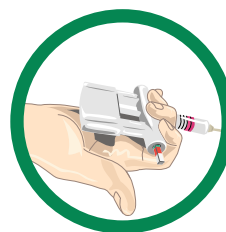

Incorrect

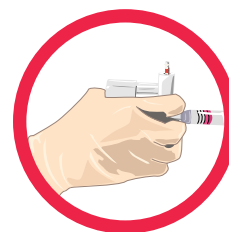

See Instructions for Use for full instructions, warnings, precautions & contraindications

Local Distributor:

clinicalinnovations.com

747 West 4170 South | Murray, Utah USA 84123
Phone +1.801.268.8200 | Fax +1.801.266.7373 |
Toll Free +1.888.268.6222

© Clinical Innovations, LLC. All rights reserved. | ART-0150 Rev.01



CLINICAL INNOVATIONS
forMOM. forBABY. forLIFE.™

*Ruler

Use ruler to measure distance from middle finger tip to interphalangeal joint and metacarpal joint.

